



Pusch®

There is no substitute for quality

TOP Grip



The TOP GRIP a product out of our instrument range with passive safety.

The intelligent construction makes the new TOP GRIP series essential instruments for the selective abrasive finishing in cases of onychosis.

TOP GRIP instruments make it easier and safer to shorten nail edges in the front and side areas or to reduce the nail plates up to the nail fold.

828T Utility model 202012009800.4

841T

805T Utility model 202012009800.4

New!

New!

New!





With our instrument range with passive safety, we offer pedicurists and podiatrists instruments with a special safety relevant construction for the practice.


A primary aim is to perform every treatment without any possible risk especially by treatments on diabetics. Our instrument range with passive safety supports the user by doing so.





By changes on infected nails the individual TOP GRIP form can shorten the overlaying nail edges safer or deformed nail plates can be selectively reduced and adapted accordingly.



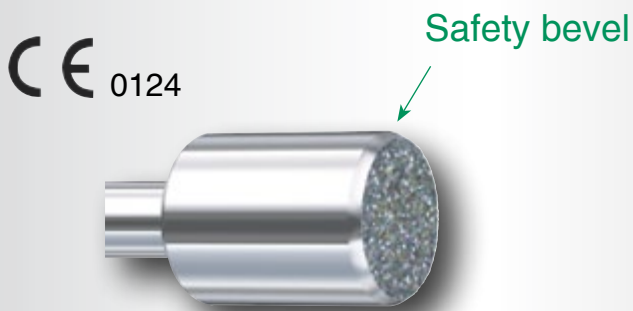

Besides the approved cylindrical basic form 840T we offer the small-sized form 841T with the option to do filigree treatments. With the new form 828T with diamond coated TOP-Surface unnaturally deformed nails and the overlaying nail in front and side areas can be reduced with a better overview and the concave safe smooth side surface fits on front and on sides of the skin contour. The new form 805T with a straight smooth surface lies supporting on the nail wall while reducing nail surface on the sides of the nails. Basically the new forms offer pedicurists/podiatrists more mobility.





Recommended maximum attainable speed: 5.000 - 10.000 rpm

Generally we recommend a decrease of rotation in the dry technique and also by treatments on high-risk patients.




			
	BUSCH	840T	8840T
	Ø ISO	050	050

			
	BUSCH	841T	8841T
	Ø ISO	032	032

			
	BUSCH	828T	8828T
	Ø ISO	065	065

			
	BUSCH	805T	8805T
	Ø ISO	065	065

 = Package unit



BUSCH & CO. GmbH & Co. KG

Postfach 1152 • 51751 Engelskirchen • Unterkaltenbach 17-27 • 51766 Engelskirchen • GERMANY
 Telefon +49 (0) 2263 860 • Telefax +49 (0) 2263 20741 • mail@busch.eu • www.busch.eu